



Ready-to-Use CO₂ Luer-Port Sensor Stick

Dear Max Mustermann


You plan to measure CO₂ in perfusion bioreactors? Online and in as little time as possible?

Then please feel invited to have a closer look at our ready-to-use CO₂ Sensor Sticks, allowing flexible integration due to a standardized luer-connector and saving time as neither set-up time nor an integration validation project are needed. An additional benefit is that the shelf-life-management of the sticks is decoupled from the disposable cultivation system. The CO₂ Sensor Sticks are particularly suitable for measurements in media with physiological osmolality.

Any questions regarding your customized application? Then do not miss to [contact our expert](#), Dr. Daniel Riechers, directly!

Your PreSens Team

Measure CO₂ in Single-use Applications

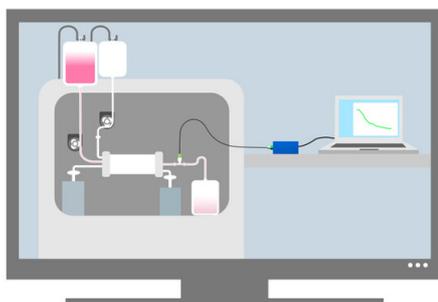


Take just 2 minutes of your time to get a detailed idea of how easy it is to integrate our small CO₂ Sensor Sticks in your customer specific flow-through cells or experimental set-up via Luer connectors.

>> [Watch our informative video](#) to see how you can benefit from our solution!

[Components needed for your set-up:](#)

- A ready-to-use **CO2 Sensor Stick** just needs to be attached to a Luer port and read-out with
- a **Polymer Optical Fiber**, which is attached to the sensor stick, allowing the non-invasive measurement. The POF is connected to
- a compact OEM solution, the Electro-Optical-Module **EOM-CO2-mini** and controlled with the
- control software **PreSens Measurement Studio PMS2** on your computer.

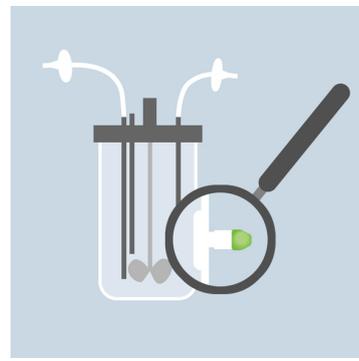

>> But please note that the SSTs should not be used in direct contact with human beings!

You would like to learn even more about PreSens Precision Sensing? Please visit our homepage www.presens.de and don't hesitate to contact us. Any feedback will be appreciated.

With kind regards

Christina Schlauderer
Communications



PreSens Precision Sensing GmbH
Am BioPark 11 - 93053 Regensburg - Germany
Phone +49 941 942 72 100, Fax +49 941 942 72 111
christina.schlauderer@presens.de, www.PreSens.de

Trade Register Ingolstadt HRB 101505, CEO: Achim Stangelmayer

[Click here to unsubscribe.](#)